

## MycoClean™ Mycoplasma Prevention Spray

# KEEP YOUR LAB SAFE AND CLEAN!

MycoClean Mycoplasma Prevention Spray is the lab's solution for decontamination.


It can be used for cell culture areas, incubators, work benches, storage containers and apparatuses.



### ● Features

- Ready-to-use / Convenient and Stable
- Non-corrosive and Non-carcinogenic
- Can be stored at RT for 2 year
- Manufactured using safe ingredients
- Elimination range: Mycoplasma, bacteria, fungi, viruses

### ● Technical data

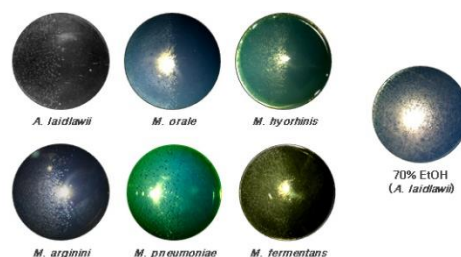

- Mycoplasma lacks a cell wall, making them unresponsive to most antibiotic treatments
- Prevention spray forms pores in cell membrane using osmotic pressure and eradicate mycoplasma.
- With this easy and convenient spray, contamination-free surfaces can be achieved with a single application.
- A test using the **6 mycoplasma species** after multiplying them on the plate and compared with 70% EtOH (For each dish, The right half of the area was sprayed out by MycoClean).
- As a result, the area **sprayed by MycoClean does not show any colony of mycoplasma**. This test also have demonstrated that 70% EtOH which is generally used for removing is not effective to mycoplasma.

### ● Product information

Product	Cat No.	Volumes
MycoClean™ Mycoplasma Prevention Spray	21083	200mL × 1 ea

### ● Related Products

Product	Cat No.	Volumes
e-Myco™ plus Mycoplasma PCR Detection Kit	25237	48 Tests
e-Myco™ VALiD Mycoplasma PCR Detection Kit	25239	48 Tests
e-Myco™ Mycoplasma PCR Detection Kit	25235	48 Tests